

ACCURATE®

Drug of Abuse Urine EZ Split Cup

Self Contained Split Specimen Drug Testing

A comprehensive urine test cup that simultaneously screens for qualitative detection of up to 15 drugs of abuse in minutes. Featuring a dual specimen chamber, the Accurate® Drug of Abuse Urine EZ Split Cup securely channels urine specimen into separate reservoirs for initial screening and retention.

Benefits

- 510(k) Cleared for Over-the-Counter use
- CLIA-Waived*
- 24-Month shelf life from date of manufacture
- Simultaneously test up to 15 drugs
- On-site drug screening
- Split specimen testing
- Qualitative results in 5 minutes
- Operator controlled test activation
- Specimen temperature indication
- Integrated specimen validity tests to detect urine adulteration

*Certain models are CLIA Waived



Features



Leak Resistant Lid

Ensure conservation of specimen with hinged, latching lid to secure urine specimen.

Large, Angled Opening

Reduce unintended mess during collection with accommodating diameter and angle of opening.



Dual Specimen Chamber

Compartmentalization of test sample and retain specimen.


Residual Sample Testing

Segregated urine specimen with minimum volume of 30 mL to a maximum of 110 mL.



Rotating Key

On demand test activation is triggered by turning unidirectional base.

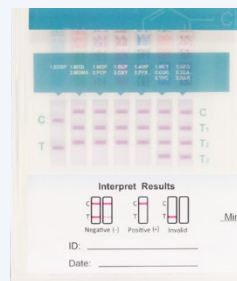

Multiplex Testing

Reduced size of result window for convenient result interpretation.

Sample-Free Result Window

Internal channels within the test cup prevent sample obscuring test results.

Flat Result Window

Results, test ID, and date unobscured by urine allowing photocopies for record keeping purposes.



ACCURATE[®]

Drug of Abuse Urine EZ Split Cup

Test Menu

FDA 510k Cleared for Over-the-Counter (OTC)		
Test		Cut-off (ng/mL)
AMP	Amphetamine	500 or 1000
BAR	Barbiturates	300
BUP	Buprenorphine	10
BZO	Benzodiazepines	300
COC	Cocaine	150 or 300
MDMA	Ecstasy	500
THC	Marijuana	50
EDDP	Methadone Metabolite	300
MTD	Methadone	300
MET	Methamphetamine	500 or 1000
MOP	Morphine	300 or 2000
OXY	Oxycodone	100
PCP	Phencyclidine	25
PPX	Propoxyphene	300
TCA	Tricyclic Antidepressants	1000

SPECIMEN VALIDITY TESTS: Creatinine (CR), Oxidants (OX), Potential Hydrogen (pH), Specific Gravity (SG), Nitrites (NI), Glutaraldehyde (GL)

**Not for screening to monitor prescription medication. Provides only a preliminary result. An alternative laboratory test must be used to confirm results. Gas chromatography/mass spectrometry (GC/MS) is the preferred method confirmation test.*

Initial Offerings

REF	Description
HCDOACI-2154S	15 Panel Split Cup with Low Cutoff BZO300, TCA1000, BAR300 MET500, COC150, THC50 AMP500, PPX300 BUP10, OXY100 OPI2000, PCP25 MTD300, MDMA500 EDDP300
HCDOACI-2155	15 Panel Split Cup BZO300, TCA1000, BAR300 MET1000, COC300, THC50 AMP1000, PPX300 BUP10, OXY100 MOP300, PCP25 MTD300, MDMA500 EDDP300

Easy to Use

