


ACCURATE FENTANYL RAPID TESTS (URINE)

A rapid, screening test for the qualitative detection of Fentanyl in human urine at a cutoff concentration of 1.0 ng/mL

3 FORMATS

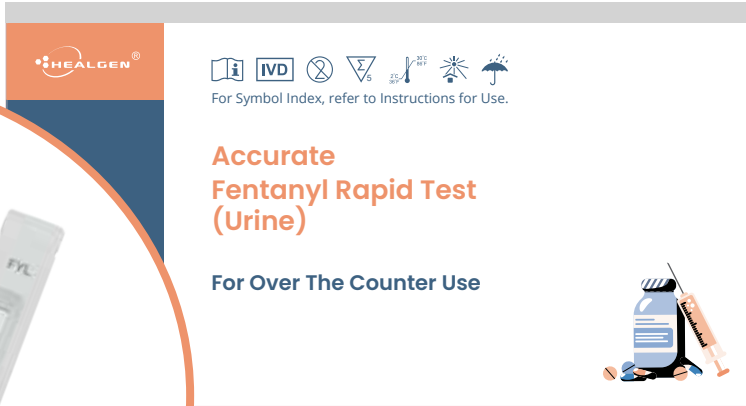

For Symbol Index, refer to Instructions for Use.

Accurate Fentanyl Rapid Test (Urine)

For Over The Counter Use

Features

- Simple collection and testing procedure
- Easy visual interpretation
- Results in 5 minutes

High Specificity for Fentanyl

Reduce False Positives: Assays that do not cross-react with norfentanyl ensure that positive results are specifically due to the presence of fentanyl, not its metabolites. This is crucial for accurately identifying current fentanyl use rather than past use.

Improve Monitoring: Specific assays help distinguish between recent fentanyl administration and the presence of residual metabolites from earlier use, aiding in more precise interpretations.


Low Cross Reactivity with the Major Urinary Metabolite Norfentanyl

*For a complete listing of this and other performance data, please see the Instructions For Use (IFU)

Fentanyl (cutoff=1.0 ng/mL)	Concentration (ng/mL)	Cross- Reactivity (%)
Fentanyl	1	100%
Norfentanyl	30,000	0.003%

Test Procedure

*This is representative of test procedures for all other formats.

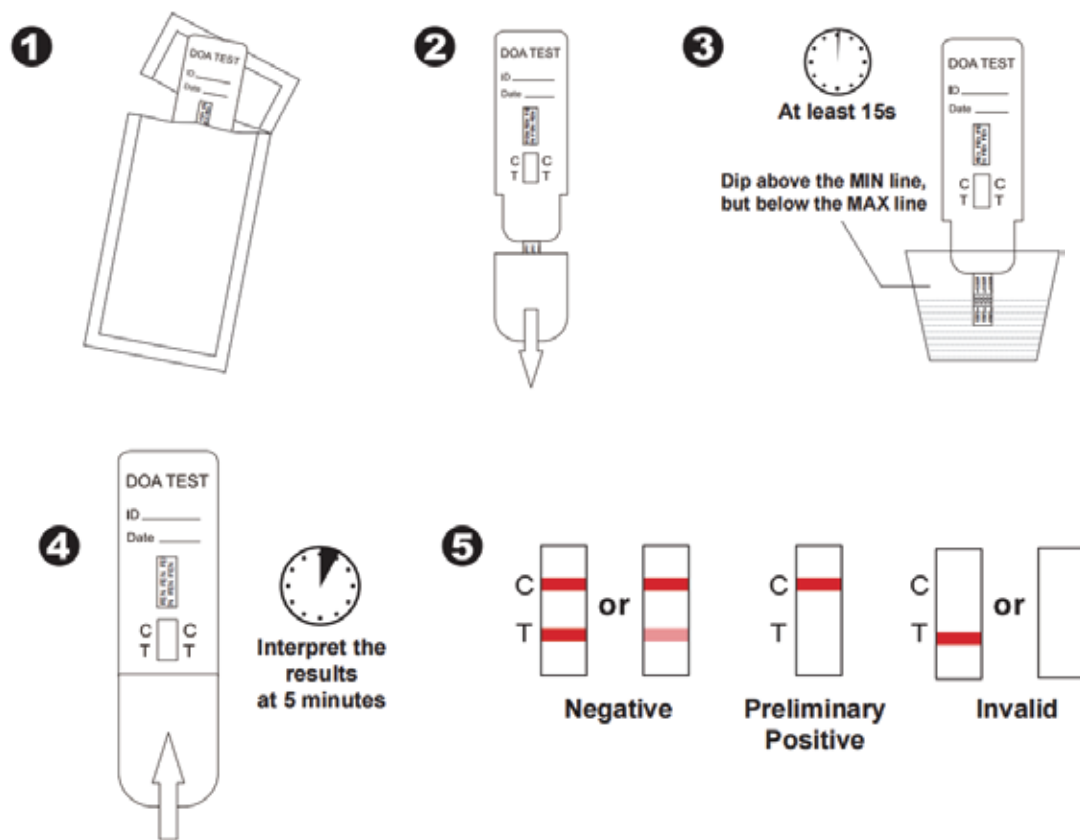

Order Information

*All tests are FDA cleared for Over-the-Counter use

Product Description	Format	Catalog No.	Kit Size
Healgen® Accurate Fentanyl Rapid Test	Strip	GBFEN-101g-H1,H2, or H5	1, 2, or 5 tests/kit
	Cassette	GBFEN-102g-H1,H2, or H5	
	Dip Card	GBFEN-105g-H1,H2, or H5	
Healgen® Accurate Rapid Fentanyl Test	Strip	GBFEN-101g	25 tests/kit
	Cassette	GBFEN-102g	
	Dip Card	GBFEN-105g	